

INFUSION PUMP

ManoPump iV1 Plus

SAFE,
PRECISE
INFUSION...
AND WARMING
IN COMPLETE
SAFETY




Practical, ergonomic and efficient

- Small and light (1.45 kg)
- Kennel mount and IV stand support included
- High-performance battery (up to 10 hours run time)
- Large touchscreen (4.3'' - 10.75 cm) for a clear display of parameters including infusion time left
- Many settings available such as:
 - ✓ Overall infusion volume
 - ✓ Requested flow
 - ✓ Infusion total time
- High sensitivity of settings (bubbles and occlusion pressure)
- Transfusion compatible



Optional fluid warmer module

Using a long fluid warmer allows the heated fluid to be brought closer to the patient, reducing heat loss, a problem frequently encountered with conventional heaters. In fact, at around ten centimetres from the outlet of the heater, there is a loss of 10 to 15°C, depending on the infusion rate and the ambient temperature.

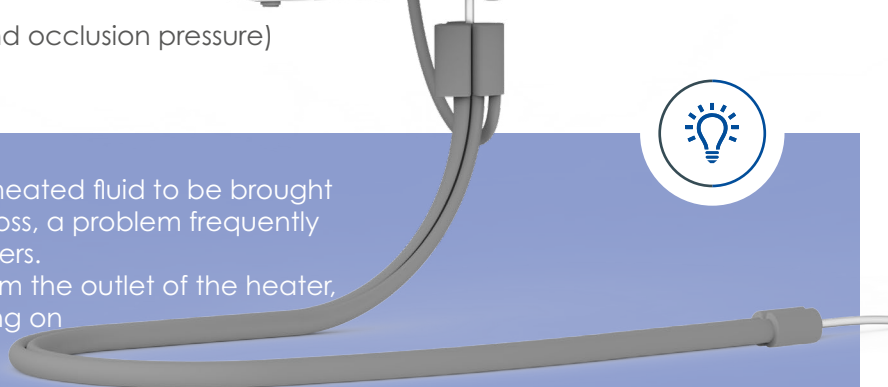

ManoPump iV1 Plus technical data:

▶ Accumulated volume	From 0.1 to 9999 ml
▶ Touchscreen	4.3'' (10.75 cm)
▶ Infusion time left	Yes
▶ Flow	From 0.1 to 2000 ml/h
▶ Modes	Flow and volume • Overall length and volume • Overall length and flow • Body mass • Saved modes • Micro mode
▶ Accuracy	+/- 5%
▶ Air bubble detection	20, 50, 100, 200, 400, 800 µl
▶ Occlusion settings	5 levels (50, 150, 300, 600, 900 mmHG)
▶ Entry capacity	Up to 2000
▶ Data review	5000
▶ Quick setting use	Access to the 20 latest settings
▶ Protection classification	IP44
▶ Dimensions and weight	13x9x14 cm - 1.45 kg
▶ Operating life	10h at 25 ml/h
▶ Power	100-240V, 50-60 Hz
▶ Safety and operating alarms	Yes
▶ Self-diagnostic	Yes
▶ Mounting system	IV stand support and kennel mount included



IV stand support



Kennel mount



22 rue Charles de Gaulle
22100 Lanvallay - France
Tel. +33 (0)2 96 85 86 76

export@manomedical.com

www.manomedical.com